

Certificate of Analysis for MRA-914G

Genomic DNA from Plasmodium falciparum, Strain FC27/PNG

Catalog No. MRA-914G

Product Description: Genomic DNA was extracted from a preparation of *Plasmodium falciparum* (*P. falciparum*), strain FC27/PNG (Papua New Guinea).

Lot¹: 63937221 Manufacturing Date: 16MAR2016

TEST	SPECIFICATIONS	RESULTS
Sequencing of Merozoite Surface Protein 2 (MSP2) gene (~ 770 base pairs)	Consistent with P. falciparum	Consistent with P. falciparum
Agarose Gel Electrophoresis	High molecular weight chromosomal DNA	High molecular weight chromosomal DNA (Figure 1)
Concentration by PicoGreen® Measurement	Report results	16.3 ng/μL
Amount	Report results	~ 500 ng per vial
Functional Activity by PCR Amplification MSP2 locus ²	~ 600-900 base pair amplicon	~ 900 base pair amplicon
OD ₂₆₀ /OD ₂₈₀ Ratio	1.6 to 2.1	2.0
Protozoan Inactivation Human erythrocytes exposed to 10% of total yield of MRA-914G³	No parasitemia observed	No parasitemia observed
Mycoplasma Contamination DNA Detection by PCR	None detected	None detected

¹MRA-914G was produced from a cell culture of BEI Resources MR-MRA-914 lot 58591230. Genomic DNA was extracted using proprietary technology.

Date: 26 JUL 2016 Signature:

BEI Resources Authentication

ATCC®, on behalf of BEI Resources, hereby represents and warrants that the material provided under this certificate has been subjected to the tests and procedures specified and that the results described, along with any other data provided in this certificate, are true and accurate to the best of ATCC®'s knowledge.

ATCC® is a trademark of the American Type Culture Collection.

You are authorized to use this product for research use only. It is not intended for human use.

E-mail: contact@beiresources.org

Tel: 800-359-7370 Fax: 703-365-2898

BEI Resources


²Primer sequences and conditions for PCR are available upon request.

³15 days in complete RPMI culture medium at 37°C in sealed flasks outgassed with blood-gas atmosphere (90% N₂, 5% CO₂, 5% O₂). Complete RPMI culture medium was changed and parasitemia checked every 2 to 4 days.



Certificate of Analysis for MRA-914G

Figure 1: Agarose Gel Electrophoresis

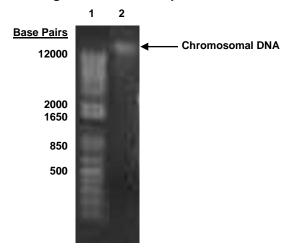

Lane 1: Invitrogen™ TrackIt™ 1 Kb Plus DNA Ladder

Lane 2: ~ 436 ng of MRA-914G

E-mail: contact@beiresources.org

Tel: 800-359-7370 Fax: 703-365-2898